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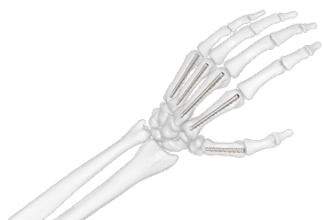
Instructions for use

ExsoMed™ INnate™ Intramedullary Threaded Nail System Hand Implant

These instructions are intended for the Operating Surgeon and supporting Healthcare Professionals. The implants and instruments are intended only for use by a licensed and trained health care professional.

DEVICE DESCRIPTION

The INnate™ Intramedullary Threaded Nail System Hand Implant is a stainless steel, fully threaded nail designed for minimally invasive repair or reconstruction of small bone fractures. It is intended for percutaneous insertion into the intramedullary cavity. The implant procedure requires three accessory devices, provided in a separate procedure kit:

- A guidewire for positioning the implant,
- A depth gauge,
- A cannulated drill to be used over the guidewire for preparing the intramedullary cavity, and
- A custom cannulated T-10 driver for implant insertion.

The INnate™ Intramedullary Threaded Nail System Hand Implant is available in 3.6mm and 4.5mm diameters, ranging from 25mm to 75mm in length with size selection based on the length of the bone being repaired.

Components

The INnate™ Intramedullary Threaded Nail System Hand Implant includes the following components to facilitate implantation:

EXINN913600, INnate Instrument Kit for 3.6mm

Components	Description	Material
10000-0303	Assembly Driver Torx T10	Assembly
10000-0002	0.045" Guide Wire Single Trocar	Stainless Steel
INN-062-4560	0.045" Guide Wire Dual Trocar	Stainless Steel
INN-031-2925	Stepped Cann. Drill Bit, 3.6mm INnate	Stainless Steel
INN-080-3575	INnate Disposable Depth Gage	Stainless Steel

EXINN914500, INnate Instrument Kit for 4.5mm

Components	Description	Material
10000-0303	Assembly Driver Torx T10	Assembly
10000-0002	0.045" Guide Wire Single Trocar	Stainless Steel
INN-062-4560	0.045" Guide Wire Dual Trocar	Stainless Steel
INN-031-3631	Stepped Cann. Drill Bit, 4.5mm INnate	Stainless Steel
INN-080-3575	INnate Disposable Depth Gage	Stainless Steel

INTENDED PURPOSE

The INnate™ system is intended for intramedullary fixation of small bone fractures, such as fractures of the metacarpals of the hand. The instrument kit is designed to facilitate placement of the implant.

INDICATIONS

The INnate™ system is indicated for fixation of intra-articular and extra-articular fractures and non-unions of small bones and small bone fragments; arthrodesis of small joints; and osteotomies of the metacarpals.

CLINICAL BENEFIT

The intended clinical benefit of the INnate™ Intramedullary Threaded Nail System Hand Implant is to provide fixation for patients with fractures as a result of trauma, osteotomy or reconstruction of bones of the appropriate size. The indirect clinical benefit of the encompassed instrumentation is to facilitate implantation protocols for fracture fixation.

TARGET POPULATION AND INTENDED USERS

The INnate™ Intramedullary Threaded Nail System Hand Implant is intended for use in skeletally mature, adult patients with intra-articular and extra-articular fractures, non-unions of small bones and small bone fragments, fractures of the hand, and osteotomies in the hand. This system is intended for use by surgeons with orthopedic training and knowledge of the indications and techniques required for fracture fixation. The device is to be implanted by the surgeon in a sterile operating room setting.

CONTRAINDICATIONS

The INnate™ system is subject to the following contraindications and are not intended for:

- Active or latent infection or marked inflammation of the treatment area.
- Sepsis.
- Insufficient quantity or quality of bone, osteoporosis, or in patients with certain metabolic diseases.
- Patients with confirmed material sensitivity.
- Patients who are unwilling or incapable of following post-operative care instructions.
- Interference or soft tissue fixation.
- These devices are not intended for screw attachment or fixation to the posterior elements (pedicles) of the cervical, thoracic, or lumbar spine.

WARNINGS and PRECAUTIONS

Warnings

- Treatment or implant may fail, including sudden failure, because of:
 - Loose fixation and/or loosening
 - Stress, including stress from inappropriate bending of the implant during surgery
 - Stress concentrations
 - Stress of weight bearing, load bearing, or excessive activity

Failure is more likely if the implant experiences increased loads due to delayed union, nonunion, or incomplete healing. Failure is more likely if the patient does not follow post-operative care instructions.

- Nerve or soft tissue damage may result from surgical trauma or the presence of an implant.
- Instrument breakage or damage, as well as tissue damage, may occur when an instrument is subjected to excessive loads, excessive speeds, dense bone, improper use, or unintended use.
- Devices of dissimilar material should not be used together in or near the implant site. Dissimilar metals in contact with each other can accelerate the corrosion process due to galvanic corrosion effects.
- Fractured implants should be removed from patient during surgery. If unable to remove, notify the patient.
- Implants may cause distortion and/or block the view of anatomic structures on radiographic images.

Precautions

- Implants and instruments are intended only for professional use by a licensed physician.
- This device requires surgical implantation under sterile conditions.
- Use devices and instruments in accordance with the Instructions for Use.
- Do not use the sterile product past the use-by date. Refer to the device label.
- INnate™ implants and accessories are single-use only. Do not use or re-sterilize an implant provided in sterile packaging if the package has been damaged. The sterility may be compromised, and the cleanliness of the implant may be uncertain. Report damaged packaging to your distributor or Acumed.
- Use of this system has not been evaluated in children or individuals that are not skeletally mature.
- Safety and effectiveness of use in patients with narrow diameter long bones has not been established.
- Inspect all components preoperatively to ensure utility. Do not attempt a surgical procedure with faulty, damaged, or suspect instruments or implants. Alternate fixation methods should be available intraoperatively.
- Reprocessing and/or reuse of single use devices may result in infection/cross contamination, and/or sudden failure due to previous stresses.
- Use of implant components from different manufacturers has not been evaluated.
- Devices may be contaminated with biohazardous materials and/or sharp materials. Observe hospital procedures, practice guidelines, and/or government regulations for the proper handling of biohazardous material and disposal of sharp materials.
- Avoid overtightening the INnate™ implant; this may strip the drive feature, preventing removal of the implant if later desired.
- When placing the guidewire, it is important to ensure it is not bent or damaged.

POTENTIAL ADVERSE EFFECTS

- Anesthesia-related problems, problems with positioning of the patient (e.g., nausea, vomiting, neurological impairments, etc.) thrombosis, embolism, infection, or injury to other critical structures such as blood vessels, excessive bleeding, etc.
- Pain, discomfort, or abnormal sensations, nerve or soft tissue damage, necrosis of bone or tissue, bone resorption, or inadequate healing from the presence of an implant or due to surgical trauma.
- Swelling, abnormal scarring, impairment in musculoskeletal function, hardware prominence, migration, loosening, bending, or breakage of the implant, malunion, nonunion or delayed union due to prolonged loading or excessive forces which may lead to implant failure and reoperation.
- Metal sensitivity, histological, allergic, or adverse foreign body reaction resulting from implantation of a foreign material. Consult our document “Metal Sensitivity Statement” at <http://www.acumed.net/ifu>.
- Injury to user.

Important: Any serious incident that has occurred in relation to the device should be reported to Acumed (through customercomplaints@acumed.net or +1.888.627.9957) and the Competent Authority of the Member State in which the user and/or patient is established.

SURGICAL TECHNIQUE GUIDE

Acumed offers one or more Surgical Technique Guides to promote the safe and effective use of this system. Consult our Surgical Technique Guides at <https://www.acumed.net/resources/documents/>. The Surgical Technique Guides also include information about the installation and calibration of the device.

Important: Surgical Technique Guides may contain important safety information.

Important: The instruments and implants in this system are intended to be used by suitably trained and qualified surgeons in a hospital operating room setting. Before treatment, the surgeon is advised to read and fully understand all instructions and communicate to the patient any relevant medical information provided therein, including the use, limitations, risks (safety communications), and possible adverse effects of the proposed treatment.

Consult the most recent versions of the Instructions for Use and Surgical Technique Guides as they are subject to change. Contact Acumed or an authorized agent for any additional information.

MRI SAFETY INFORMATION

Many Acumed implants have been evaluated for safety in the MR environment. Consult our publication “Acumed Implants in the MR Environment” at <http://www.acumed.net/ifu>.

LIFETIME

- Once installed, implants are expected to provide fixation, fusion and physiological support and have an effective life during bone healing. The implants are biocompatible and may remain implanted at the discretion of the surgeon or patient.

STERILITY

- The INnate™ Intramedullary Threaded Nail System Hand Implant and accessories are provided **STERILE** via gamma irradiation, in accordance with ISO 11137, to achieve a Sterility Assurance Level (SAL) of 10^{-6} .
- **DO NOT USE IF THE STERILE PACKAGE IS DAMAGED. DO NOT USE AFTER EXPIRATION DATE.**

IMPLANT

MATERIALS

- The INnate™ implant is manufactured from 316L stainless steel for use in surgical implants per ASTM F138.
- Consult our document “Metal Sensitivity Statement” at <http://www.acumed.net/ifu> for a detailed description of the material composition of Acumed metal implants.

SINGLE USE IMPLANTS

- Implants are intended for single use only, as indicated on the label. This product has not been validated for multiple uses. Reusing single use implants may increase the risk of failure and cross- contamination.
- Dispose of any unused implant that is contaminated with human blood or tissue in accordance with hospital procedures, practice guidelines, and/or government regulations for the proper handling of biohazardous material. Do not process a contaminated implant.

IMPORTANT

- Physiological dimensions limit implant sizes. Select the type and size of implant that best meets the patient’s requirements for close adaptation and firm seating with adequate support.
- Implants are designed to hold fractured bones in place during bone healing at the site of fracture. Implants are not designed for excessive loadbearing. Improper selection or improper implantation of the device may increase the possibility of loosening or migration.
- Do not bend the implant except as indicated in the corresponding Surgical Technique Guide. Repeated or excessive bending may weaken the implant and cause failure at a later time.
- Only combine implants when they are intended for that purpose.
- Protect implants against scratching and nicking to prevent stress concentrations, which can result in failure.
- Prevent unused implants from becoming soiled.
- The color of anodized implants may change over time due to processing. This color change does not affect the mechanical properties of the implants.

INSTRUMENTS

MATERIALS

The instruments are manufactured from various grades of stainless steel and anodized aluminum.

MULTIPLE USE and SINGLE USE

The INnate™ instruments are intended for single use only.

Single Use Instruments

- Do not reuse single use instruments. Reuse or reprocessing of devices labeled as “single use” can result in decreased mechanical and clinical performance of the devices which may result in patient harm.
- Reuse or reprocessing of single use instruments may create a risk of contamination (e.g., due to the transmission of infection material from one patient to another), which may result in patient harm.
- Dispose of single use instruments after use on a single patient, during a single procedure.

IMPORTANT

- Protect instruments against scratching and nicking to prevent stress concentrations, which can lead to instrument failure.

STORAGE CONDITIONS

STORAGE OF PACKAGED STERILE PRODUCT

Devices should be stored in an area that provides protection from dust, pests, and temperature/humidity extremes.

SAFE DISPOSAL

Devices may be contaminated with biohazardous materials and/or sharp materials. Observe hospital procedures, practice guidelines, and/or government regulations for the proper handling of biohazardous material and disposal of sharp materials.

SUMMARY OF SAFETY AND CLINICAL PERFORMANCE

The Summary of Safety and Clinical Performance Report for the implants may be obtained from the European dataset on medical devices (EUDAMED) at <https://ec.europa.eu/tools/eudamed>. Upon upload to EUDAMED, the search terms to view this document will be:

- INnate™ implants: 0806378BUDI71Q9
- INnate™ 3.6 instruments: 0806378BUDI72QF
- INnate™ 4.5 instruments: 0806378BUDI72QH

Important: Any serious incident that has occurred in relation to the device should be reported to Acumed and the competent authority of the Member State in which the user and/or patient is established. To report a complaint about this product, call Customer Service at +1.888.627.9957 or send an email to customercomplaints@acumed.net.

Symbols Glossary

Symbol	Description	ISO 15223-1
 www.acumed.net/ifu	Consult the electronic instructions for use (eIFU) at www.acumed.net/ifu	5.4.3
STERILE R	Sterilized using irradiation	5.2.4
	Single sterile barrier system with protective packaging inside	5.2.13
	Use-by date	5.1.4
REF	Catalogue number	5.1.6
LOT	Batch code	5.1.5
EC REP	Authorized representative in the European Community / European Union	5.1.2
MD	Medical device	5.7.7
	Manufacturer	5.1.1
	Date of manufacture	5.1.3
	Do not re-sterilize	5.2.6
	Do not re-use	5.4.2
UDI	Unique Device Identifier: Indicates a carrier that contains Unique Device Identifier information	5.7.10
	Do not use if package is damaged and consult instructions for use / do not use if the product sterile barrier system or its packaging is compromised.	5.2.8
	Non-pyrogenic	5.6.3
	Keep Dry: Protect from moisture. Keep away from rain.	5.4.3
	MR Conditional: An item with demonstrated safety in the MR environment within defined conditions. (ASTM F2503 7.3.2)	
	The reticle is a registered trademark of Acumed. It may appear alone or with the Acumed name.	
	CE marking of conformity, Article 17 of EU Directive 93/42/EEC or Article 20 of Regulation (EU) 2017/745. CE marking may be accompanied by the identification number of the notified body responsible for conformity assessment.	

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Implante para mano del sistema de clavos roscados intramedulares ExsoMed™ INnate™

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Fax: +1.503.520.9618

www.acumed.net

These materials contain information about products that may or may not be available in any particular country or may be available under different trademarks in different countries. The products may be approved or cleared by governmental regulatory organizations for sale or use with different indications or restrictions in different countries. Products may not be approved for use in all countries. Nothing contained in these materials should be construed as a promotion or solicitation for any product or for the use of any product in a particular way that is not authorized under the laws and regulations of the country where the reader is located. Nothing in these materials should be construed as a representation or warranty as to the efficacy or quality of any product, nor the appropriateness of any product to treat any specific condition. Physicians may direct questions about the availability and use of the products described in these materials to their authorized Acumed distributor. Specific questions patients may have about the use of the products described in these materials or the appropriateness for their own conditions should be directed to their own physician.

PKGI-014516 Rev C | Effective 1-2026 | © 2023 Acumed® LLC

ExsoMed™ and INnate™ are Trademarks of ExsoMed Corporation. ExsoMed™ Corporation is a wholly owned subsidiary of Acumed LLC. Acumed® is a registered trademark of Acumed LLC.

Subsidiaries:

Acumed Ltd

Huebner House
The Fairground
Andover
Hampshire UK SP11 0QN
Tel: +44 1264 774450

Acumed Iberica

C/ Proción 1
Edificio Oficor
28023 Madrid – Spain
Tel: +34 913516357

Acumed Beijing

Room A1206, Horizon International
Tower
No. 6, Zhichun Road
Haidian District
100088 Beijing, China
Tel: +86 10 82001303

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Acumed GmbH