


# EXTREMIFIX™

Cannulated Screw System | Mini & Small

## Product Information

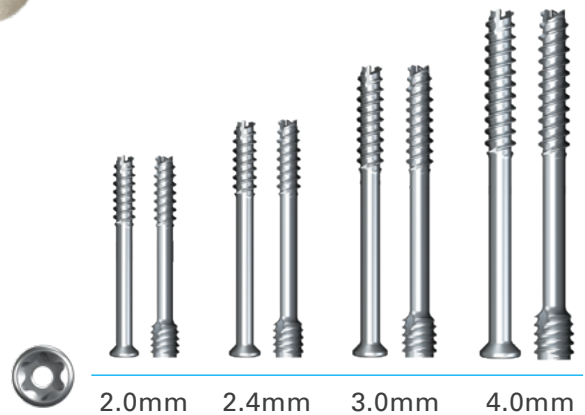

*Rethinking Possibilities, Reshaping Lives*

# ExtremiFix Mini & Small

The ExtremiFix Mini & Small Cannulated Screw System offers multiple screw and instrument options and sizes for fracture fixation, fusion, and osteotomy procedures of the upper and lower extremities. The system offering includes Ø2.0mm, Ø2.4mm, Ø3.0mm, and Ø4.0mm partially-threaded, headed and headless screws in various lengths.

## Headed Screw Design

## Headless Screw Design

**Low-Profile Head**  
Minimizes soft tissue irritation

**Conical Head**  
Minimizes hardware prominence after countersinking and creates concentric compression forces



**Hexalobe Drive Recess T7 and T10**  
Higher torque resistance



**Reverse Cutting Flutes**  
Easy removal

**Double Lead**  
(Ø3.0mm & Ø4.0mm only)  
Faster advancement

**Sharp Threads**  
Bony engagement,  
purchase, stability, fixation

**Self-Drilling and Self-Tapping**  
Easy insertion

**Headless Screw Compression**  
Differential pitch between the cancellous, greater-pitched, leading threads and the cortical, finer-pitched, trailing threads

**Minimize Protrusion**  
Allows for flush placement within the bone



## Washers

The ExtremiFix Mini & Small Cannulated Screw System includes washers for use with the Ø2.0mm, Ø2.4mm, Ø3.0mm, and Ø4.0mm headed screws. Washers may be used to provide additional contact area in patients with soft bone and are laser-marked to indicate which side must be in contact with the screw head.



Refer to the provided instructions for use for the complete indications, contraindications, warnings, and instructions for use.



3885 Arapaho Rd.  
Addison, TX 75001  
Customer Service: 800.456.7779  
Outside the U.S.: 001.972.677.4600  
Fax: 800.390.2620  
Fax Outside the U.S.: 001.972.677.4709  
E-mail: [customer.service@osteomed.com](mailto:customer.service@osteomed.com)  
[www.osteomed.com](http://www.osteomed.com)