




acumed®

Stableloc External Fixation System

Product Overview

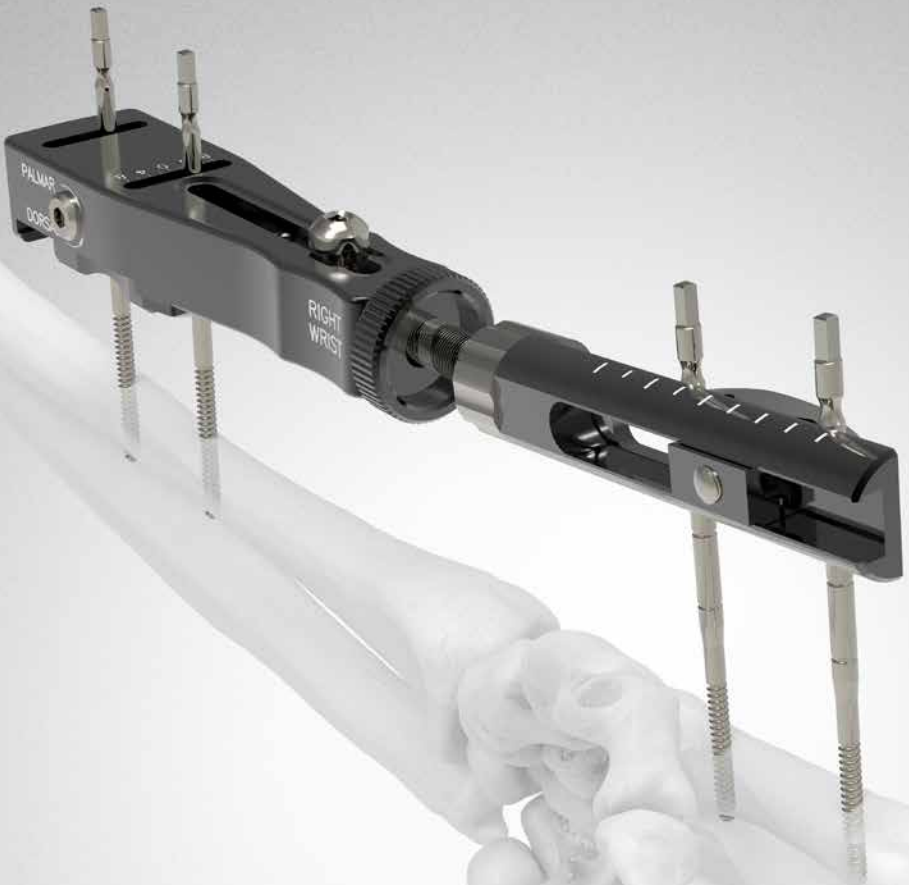

A COLSON ASSOCIATE

Acumed offers the Stableloc External Fixation System for the treatment of complex distal radius fractures. The Stableloc External Fixator is designed to hold pins in place for fracture reduction or distraction to align the wrist and to provide ligamentotaxis to the wrist on a temporary basis while the distal radius heals.


FEATURES OF THE SYSTEM INCLUDE:

Versatility allows surgeon to gain initial reduction with the fixator in place and then separately adjust only those planes needing correction. Adjustments may be made independently for distraction, radial/ulnar deviation, flexion/extension and dorsal/palmar translation. Radiolucent body aids in surgeon visualization during fracture reduction and distraction. The Stableloc External Fixator is designed to be lightweight. The sterile system includes all instruments and pins needed to complete a case and adjust the Stableloc External Fixator. 3.2 mm pins are designed as self drilling and come standard in the Stableloc kit.

The Stableloc External Fixation System may also be used with the following Acumed products, as applicable:

- Acu-Loc® Volar Distal Radius System
- Acu-Loc® 2 Wrist Plating System
- Acutrak® Headless Compression Screws
- Acutrak 2® Headless Compression Screws
- ARC Wrist Tower
- Hand Fracture System
- Osteotomy System
- SLIC Screw® System
- Small Bone External Fixation System

For more information, please visit: go.acumed.net/stableloc

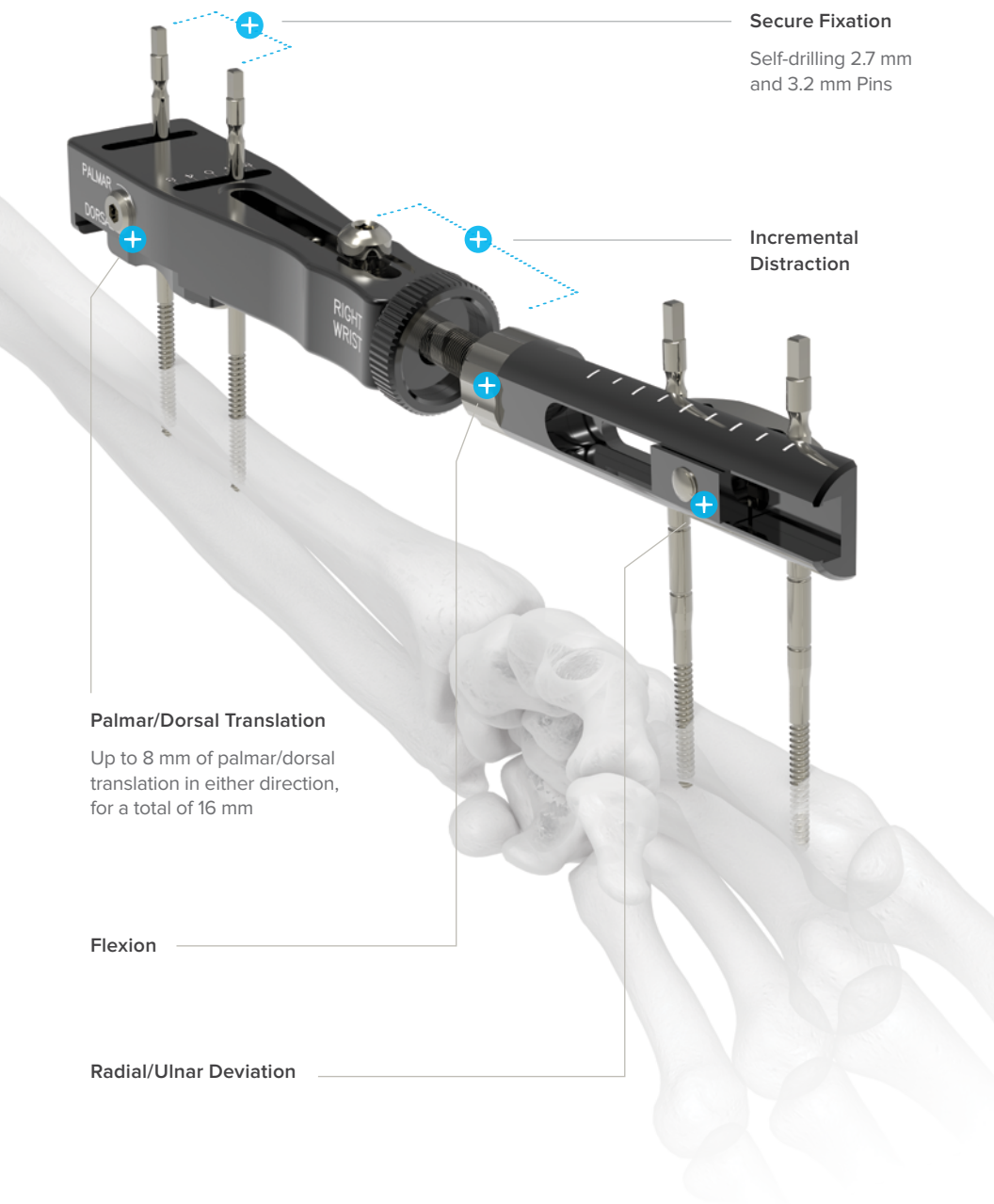

Secure Fixation

Self-drilling 2.7 mm
and 3.2 mm Pins

**Incremental
Distraction**

Palmar/Dorsal Translation

Up to 8 mm of palmar/dorsal
translation in either direction,
for a total of 16 mm

Flexion

Radial/Ulnar Deviation



Acumed Headquarters
5885 NW Cornelius Pass Road
Hillsboro, OR 97124
Office: +1.888.627.9957
Office: +1.503.627.9957
Fax: +1.503.520.9618
www.acumed.net

HNW40-12-A | Effective: 09/2015 | © 2015 Acumed® LLC

These materials contain information about products that may or may not be available in any particular country or may be available under different trademarks in different countries. The products may be approved or cleared by governmental regulatory organizations for sale or use with different indications or restrictions in different countries. Products may not be approved for use in all countries. Nothing contained on these materials should be construed as a promotion or solicitation for any product or for the use of any product in a particular way which is not authorized under the laws and regulations of the country where the reader is located. Specific questions physicians may have about the availability and use of the products described on these materials should be directed to their particular local sales representative. Specific questions patients may have about the use of the products described in these materials or the appropriateness for their own conditions should be directed to their own physician.

Acumed®, Acu-Loc®, Acutrak® and Acutrak 2® are registered trademarks of Acumed, LLC