

The SOLUTION for your THORACIC STABILIZATION needs!

RibLoc®

Rib Fracture Plating System



RibLoc® U+

Chest Wall Plating System




AcuTie® II

The Next Generation
Sternum Closure System




BioBridge®

Resorbable Chest Wall
Stabilization Plate

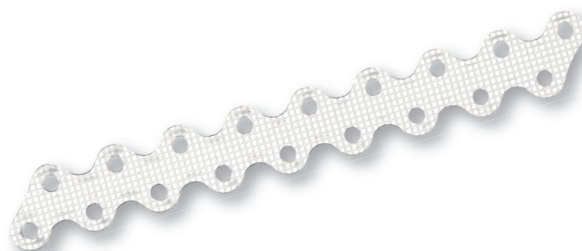

Back to health. Back to work. Back to **life.**

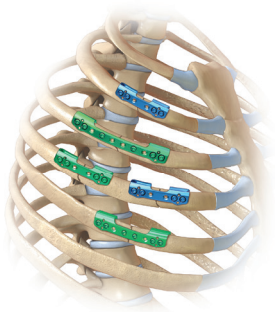

RibLoc® Rib Fracture Plating System

- MINIMALLY - INVASIVE** Plates available in three lengths with the shortest being only 46mm
- STABLE FIXATION** U-shape supports the rib on three surfaces
Avoids neurovascular bundle on inferior margin
Screws lock posteriorly
- STRAIGHT - FORWARD** Color-coded system




Ribloc® U Plus Chest Wall Plating System

- STABLE** U-clip fixation doesn't rely on bone quality
Anterior and Posterior Dual-Locking technology
- VERSATILE** Addresses diverse fracture patterns
Targets a broader range of patient needs
- CUSTOMIZABLE** U-clips compress to match rib thickness
Bending instrumentation
- EFFICIENT** Multi-purpose instrumentation to eliminate steps
One screw, drill and driver
- U Plus 90** Low profile instrumentation
90° angled handpiece with LED light
Power to compress U-clip, drill and install screws



BioBridge® Resorbable Chest Wall Stabilization Plate

- STABILIZES** Fractures or osteotomies following thoracotomy procedures
Costochondral junction fractures
Traumatic chest wall injuries
Reconstructive procedures of the thorax
- NON - PERMANENT** Polylactide acid blend offers strength
Maintains stability for 6 months
Metabolizes within 18-24 months through hydrolysis

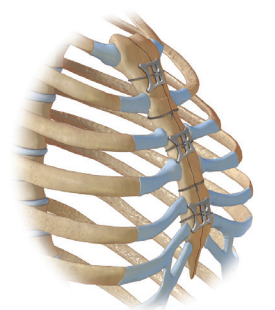

AcuTie® II Sternum Closure System

- STABILITY** Resists rostral-caudal and anterior-posterior motion
Provides approximation and compression at the osteotomy
73% increase in force required to tear through, compared to wire alone
- QUICK INSTALLATION** No measuring
No drilling or driving screws
No plate bending or cutting
No extra dissection
- COST - EFFECTIVE** Low implant cost
Reduced OR time



21421 NW Jacobson Road
Suite 700
Hillsboro, OR 97124

866.623.4137

www.acuteinnovations.com

