







The **Logic Mandibular Distraction System** is an intraoral bone distractor. It features various curved and straight bars activated with a threaded wire that has screw holes that are fixed to bone via 1.6mm or 2.0mm bone screws. The distractor is available in right and left versions. The threaded wire is activated by a hex driver and is capable of distraction lengths of up to 25mm.

- Low profile intra-oral device for use in cases where distraction osteogenesis is indicated
- Choice of three devices with fixed curvilinear trajectories that mimic the natural growth curve of the mandible
- Utilizes a patented flexible wire that provides access through soft tissue for ease of activation
- The activation wire can be cut during the consolidation period reducing the risk of infection while promoting soft tissue healing


Logic[™] **Distractors**

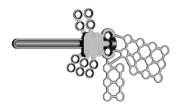

216-0119 LOGIC™ Distractor, Straight Right

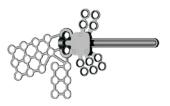


216-0118 LOGIC™ Distractor, Straight Left

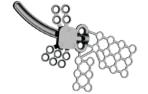

216-0113 LOGIC™ Distractor, 36mm Right



216-0112 LOGIC™ Distractor, 36mm Left

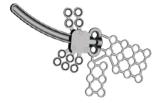

216-0111 LOGIC™ Distractor, 52mm Left



216-0110 LOGIC™ Distractor, 52mm Right

1.6mm Screws





| 208-1603 | 1.6mm x 3mm Screw |
|----------|-------------------|
| 208-1604 | 1.6mm x 4mm Screw |
| 208-1605 | 1.6mm x 5mm Screw |
| 208-1606 | 1.6mm x 6mm Screw |
| 208-1608 | 1.6mm x 8mm Screw |

AutoDrive® Screws





| 211-1604 | 1.6mm x 4mm Auto-Drive® | Screw |
|----------|-------------------------|-------|
| 211-1605 | 1.6mm x 5mm Auto-Drive® | Screv |
| 211-1606 | 1.6mm x 6mm Auto-Drive® | Screw |

OsteoMed.com 2


Activation Wires



Instruments



OsteoMed.com

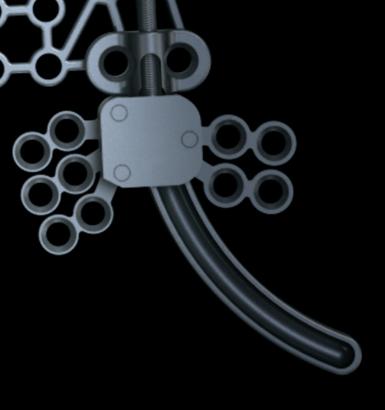



Rethinking Possibilities, Reshaping Lives

OSTEOMED

3885 Arapaho Rd. Addison, TX 75001

Customer Service: 800.456.7779 Customer Service Fax: 800.390.2620

www.osteomed.com P/N 030-1218 Rev.J



